

Health & Beauty Creator

for your healthy life,
for the wonderful world

Before & After



• Before • After • Before • After • Before • After

Technical Specifications

- Max. power consumption : 70W
- Output Power : 0~3J/cm²
- Output Frequency : 4/7 MHz
- Cartridge type : 1.5/3.0/4.5mm
- Spacing : 1.0~ 2.0 mm
- Length : Half, Full
- Lifespan : 4MHz 16,000/ 7MHz 10,000 Lines
- Power Input : AC100/240V, 50/60Hz
- Dimensions : 35X34X15cm
- Weight : 5Kg(main unit),15Kg(Trolley)

High Intensity
Focused Ultrasound

ce||thion



daeOyang
Medical

Donghwa Medical Instrument Complex, 1656-1, Donghwa-ri, Munmak-eup, Wonju-si, Gangwon-do, 220-801 Korea
Tel +82 33 747 4465~6 Fax +82 33 747 1133 E-mail daeyang@daeyangmed.com www.daeyangmed.com

daeOyang
Medical

Health &
Beauty Creator




ce||thion
High Intensity Focused Ultrasound

● What is Cellthion?

The SAFEST Ultrasound treatment system applying HIFU (High Intensity Focused Ultrasound)


- Transducers micro-focus ultrasound energy below the skin to non-invasively denature collagen.
- Tighten and actually lift skin tissue without disrupting the surface of skin by HIFU

● Cellthion Effects

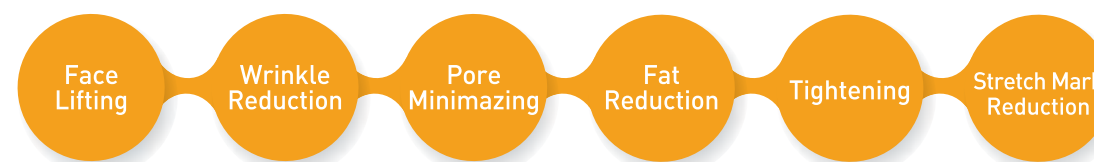

● Strengths

- Max 16,000lines
- Smart motor control
- Preventing damage of the film by double covering
- Keep air-pressure balance
- Full magnetic
- Installed Detect sensor
- Ergonomic design
- Vertical Direction Usage
- Pre-set clinical guidelines

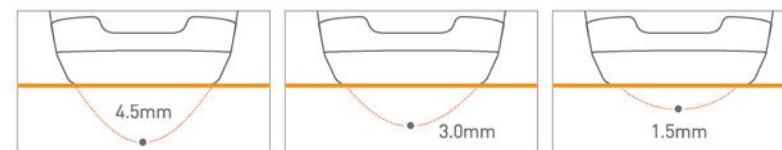

● Benefits

- Immediate Result
- Non-invasive treatment
- No downtime
- Proven ultrasound technology
- Natural, gradual results over 2-3months
- Convenient, 30minute procedure
- No special recuperation needs
- One Session

● How to Process?

- Focused ultrasound is delivered non-invasively and heats deep tissue without affecting the surface of skin.
- The heat creates triggers that is the body's natural response mechanism, which is to rejuvenate tried collagen and supplement it with fresh, new collagen.
- After that, weak collagen is reorganized and strengthened by the growth of new collagen, which tightens and lifts the skin. That is "Lifting Effect".



• Handpiece

Cartridge •